Manual de Uso

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Operating Instructions

Please read and follow the instructions in this manual carefully before use.


AlleviAider™

How to Use

1. Make sure the polarities on the 9V battery and the battery chamber match up. Open the battery cover and insert the 9V battery into the battery chamber. Close the battery cover by sliding it back to the original position.
2. Make sure to install the device on the adjustable knob before turning on the power. Point the aperture at the advised location (described in this manual).
3. To avoid direct eye contact with the light source, please put on the safety goggles provided before use.
4. Each session is 10 minutes. The device will power off automatically and make a "beeping" sound at the end of the session.
5. The optimal distance between the device and the skin is about 1 centimeter.

Caution

1. To avoid direct eye contact with the light, please put on the safety goggles provided before use.
2. If you experience any discomfort while using this product, turn the device off immediately and contact your healthcare provider.
3. In case of emergency, press the "Power ON/OFF" switch again to stop the treatment.
4. Keep this product out of children's reach to avoid any risk.
5. Avoid accidental shocking or impact to this product.
6. If you won't be using this product for a long period of time, please remove the battery and storage the device in its original box to prevent damage.
7. Do not use this device in high temperatures, the bathroom or other locations with high humidity.
8. Do not use this device while driving or sleeping.
9. Do not use this device on infants or people incapable of expressing their feelings.

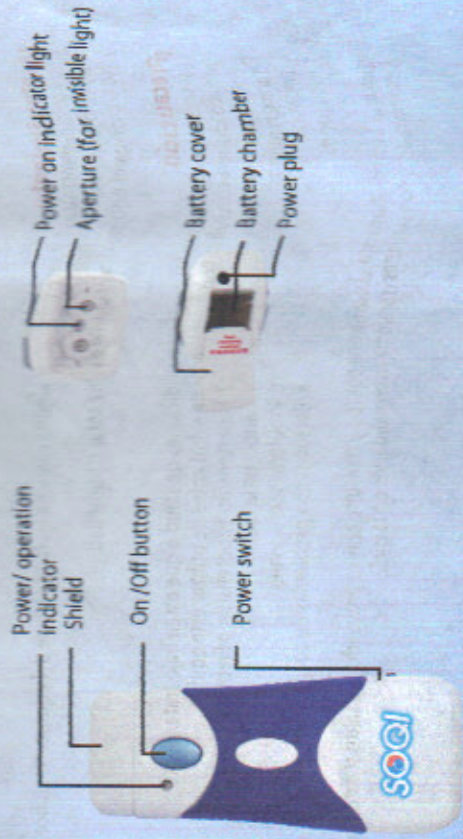Usage and Precautions

The effect is cumulative. The effect is cumulative. The effect is cumulative.

Warning

- To prevent fire or electric shock, do not use the device in the following situations:
 1. In a wet area.
 2. In a high humidity area.
 3. In a high temperature area.

Maximum Use

1. Use the device for no more than 10 minutes at a time.
2. Avoid repeated use of the device.

Caution

1. Avoid direct eye contact with the light source.
2. Avoid direct eye contact with the light source.
3. Avoid direct eye contact with the light source.
4. Avoid direct eye contact with the light source.
5. Avoid direct eye contact with the light source.
6. Avoid direct eye contact with the light source.
7. Avoid direct eye contact with the light source.

*Immediately stop using the device if you experience any discomfort or pain.

Operating Instructions

Please read and follow the instructions in this manual carefully before use.


AlleviAider™

How to Use

1. Make sure the polarities on the 9V battery and the battery chamber match up. Open the battery cover and insert the 9V battery into the battery chamber. Close the battery cover by sliding it back to the original position.
2. Make sure to install the device on the adjustable belt before turning on the power. Point the aperture at the advised location (described in this manual).
3. To avoid direct eye contact with the light source, please put on the safety goggles provided before use.
4. Each session is 10 minutes. The device will power off automatically and make a "beeping" sound at the end of the session.
5. The optimal distance between the device and the skin is about 1 centimeter.

Caution

1. To avoid direct eye contact with the light, please put on the safety goggles provided before use.
2. If you experience any discomfort while using this product, turn the device off immediately and contact your health-care provider.
3. In case of emergency, press the "Power ON/OFF" switch* again to stop the treatment.
4. Keep this product out of children's reach to avoid any risks.
5. Avoid accidental shocking or impact to this product.
6. If you won't be using this product for a long period of time, please remove the battery and storage the device in its original box to prevent damage.
7. Do not use this device in high temperatures, the bathroom or other locations with high humidity.
8. Do not use this device while driving or sleeping.
9. Do not use this device on infants or people incapable of expressing their feelings.

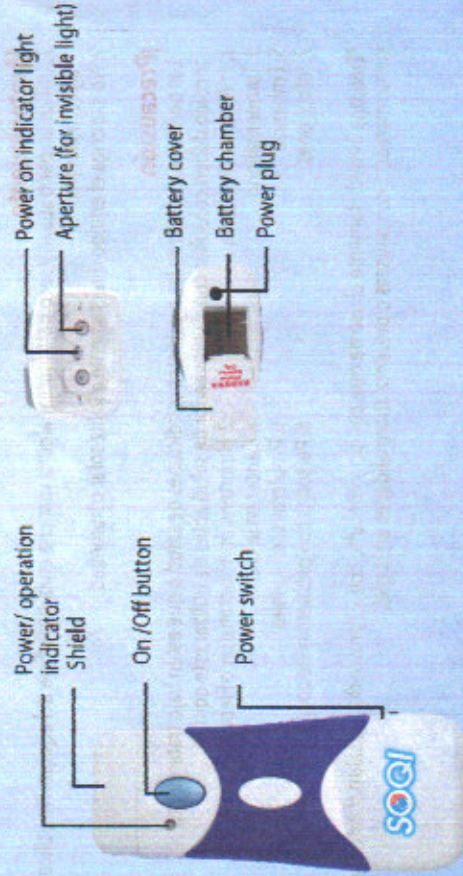Wavelength 830nm ±10nm

Irradiation Range 186 mm²

Total Output Power 2 x 30mW

Power Input: 100-240V AC (or)
9A Alkaline Battery
Output: 9V DC, 1A

Output impulse frequency 10Hz

Automatic Power OFF Approximately 10 minutes

Dimension 4.72inch x 2.047inch x 1.10inch

Indicators
Red: Emitting direction of light
Green: Starting and preparing
Orange: Low power/battery need to be replaced

Accessories Adjustable belt, Goggles

Usage and performance

This product uses the energy emitting from infrared radiation for the purpose of alleviating pain, which facilitates local blood circulation, eases nerve, muscle, shoulder and/or arthritis pain.

Warning

*To prevent malfunction, do not use this device with the following electronic medical instruments:

1. Implanted device such as cardiac pacemaker
2. Life-sustaining device such as cardiovascular/pulmonary pacemaker.
3. Wearable device, such as ECG machine.

*Accidents might happen when used with the above mentioned devices at the same time.

*Stop using this product immediately if it starts to function incorrectly and contact Hsin Ten Enterprise U.S.A for repairs.

*Do not repair, disassemble or modify this device by yourself. Please contact Hsin Ten Enterprise U.S.A for repairs.

*Do not insert any foreign object into this device.

Maintenance

1. Use dry, soft cloth to wipe any stain on the surface of the device. Do not use water or alcohol as a cleaner.
2. Avoid exposing this product directly under the sun or in humidity.

Caution

*Patients with the following conditions or under medical treatment should consult their health care provider before using this device.

1. Acute disease
2. Hipertension/ Hypotension or fever
3. Malignant tumor
4. Skin Cancer
5. During pregnancy
6. Sensory abnormalities
7. Infectious diseases
8. People with coagulation disorders.

*Immediately stop the usage of this product and contact your health care provider if you feel dizzy, nauseated or present any skin abnormalities.